



# Microgravity Biomaterials

---



## Overview of the NASA Microgravity Biomaterials Program

by

James Patton Downey  
EM31/MSFC



# Microgravity Biomaterials Sponsor

---



Microgravity Biomaterials Sponsor: The microgravity research program is sponsored by the Space Life and Physical Sciences Division of the Human Exploration and Operations Directorate at NASA Headquarters

- Director: Marshall Porterfield
- Program Executive for Physical Sciences: Fran Chiaramonte
- Biomaterials research project management is delegated to MSFC

Six areas of microgravity research are supported

- Biophysics
- Combustion
- Complex Fluids
- Fluids
- Fundamental Physics
- Materials



# Reasons for Performing Microgravity Studies

---



Microgravity is a unique environment for laboratory studies

- A lab experiment on the International Space Station is in a state of nearly perfect ballistic freefall
- Note: slight atmospheric resistance and the offset of the location of the attachment of the experiment from the ISS center of mass prevent a perfect ballistic freefall




# Reasons for Performing Microgravity Studies

---



Near elimination of buoyancy driven convection

- e.g. Diffusion dominated growth in directional Solidification experiments on crystalline and dendrite formation – Studies by Dr. Trivedi of Iowa State University

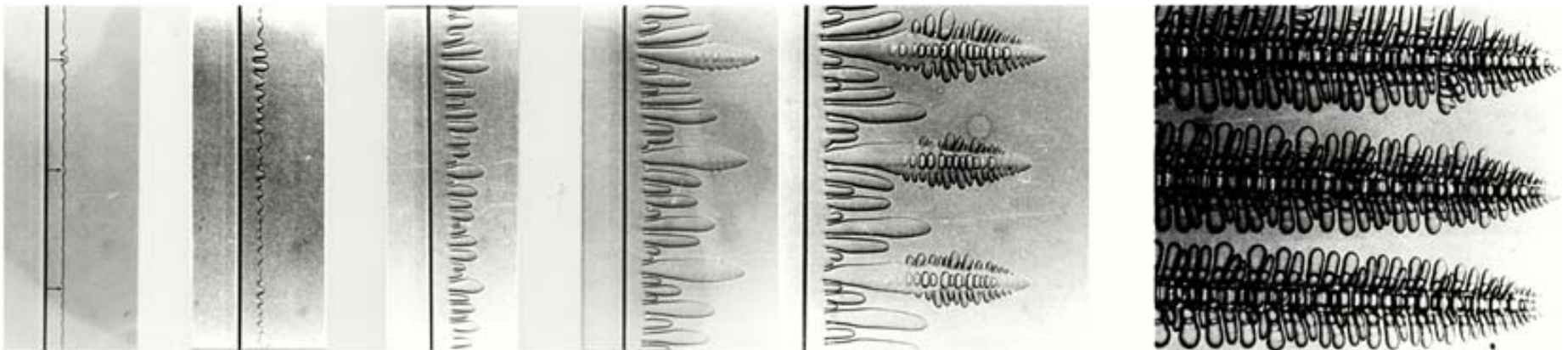



# Reasons for Performing Microgravity Studies

---



Near elimination of pressure head effects

- e.g. Stress-strain studies of granular samples at low-confining pressures – The Mechanics of Granular Materials investigations by Stein Sture of the University of Colorado






# Reasons for Performing Microgravity Studies

---



Near elimination of sedimentation

- e.g. Coarsening of spherical droplets (molten Sn in Pb) via Ostwald Ripening – Microgravity Studies by Dr. Peter Voorhees of Northwestern University

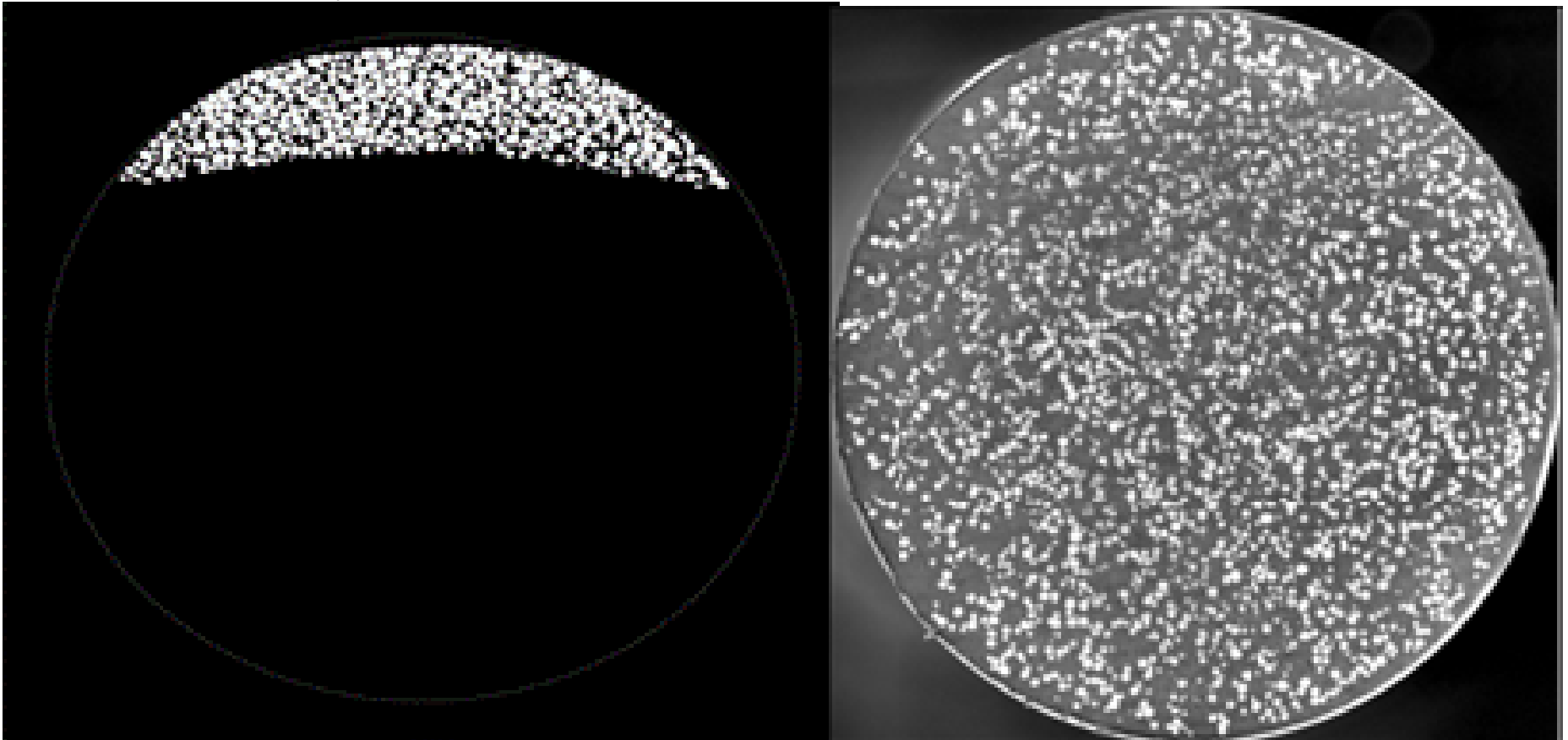



# Reasons for Performing Microgravity Studies

---



Note: Studies of molecular level interactions are typically not a reason for performing microgravity experiments

- Gravity is so weak that electrostatic interactions between molecules overwhelm any gravitational effects



# Microgravity Biomaterials Science Selection

---



## SLPS Science Selection Process:

Traditionally projects are selected from proposals submitted in response to a NASA Research Announcement soliciting response in a defined area of interest

- An NRA soliciting proposals to “mine” data from previous spaceflight experiments is planned in this year.
- Data from previous microgravity physical sciences experiments sponsored by SLPS and conducted on the ISS is being collected for a database accessible to the user community
  - POC: [Ben.Henrie@NASA.Gov](mailto:Ben.Henrie@NASA.Gov), 256-544-2446.

Alternatively scientific project requirements may be determined by a Science Definition Team with post-flight analysis of results conducted by a research group to be determined near or after the flight.





# Microgravity Biomaterials History

---



The precursor to the current biomaterials program was a series of experiments in which protein crystals, crystals of viruses, and nucleic acid crystals were grown in microgravity conditions by a variety of methods


The investigations were largely conducted by Drs. Alexander McPherson, Larry Delucas, and Daniel Carter with associated colleagues

Improved crystalline order from a number of samples such as insulin complexes, satellite tobacco mosaic virus, etc. allowed more detailed structural data from diffraction patterns to be obtained.

In a number of cases crystals large enough for neutron diffraction studies were obtained.

Reasons for the improvement attributed to factors such as slower growth, coarsening/annealing, and/or improved partitioning of large impurities from the growing crystal surface.


# New and Improved Published Macromolecular Structures Resulting from Microgravity Research

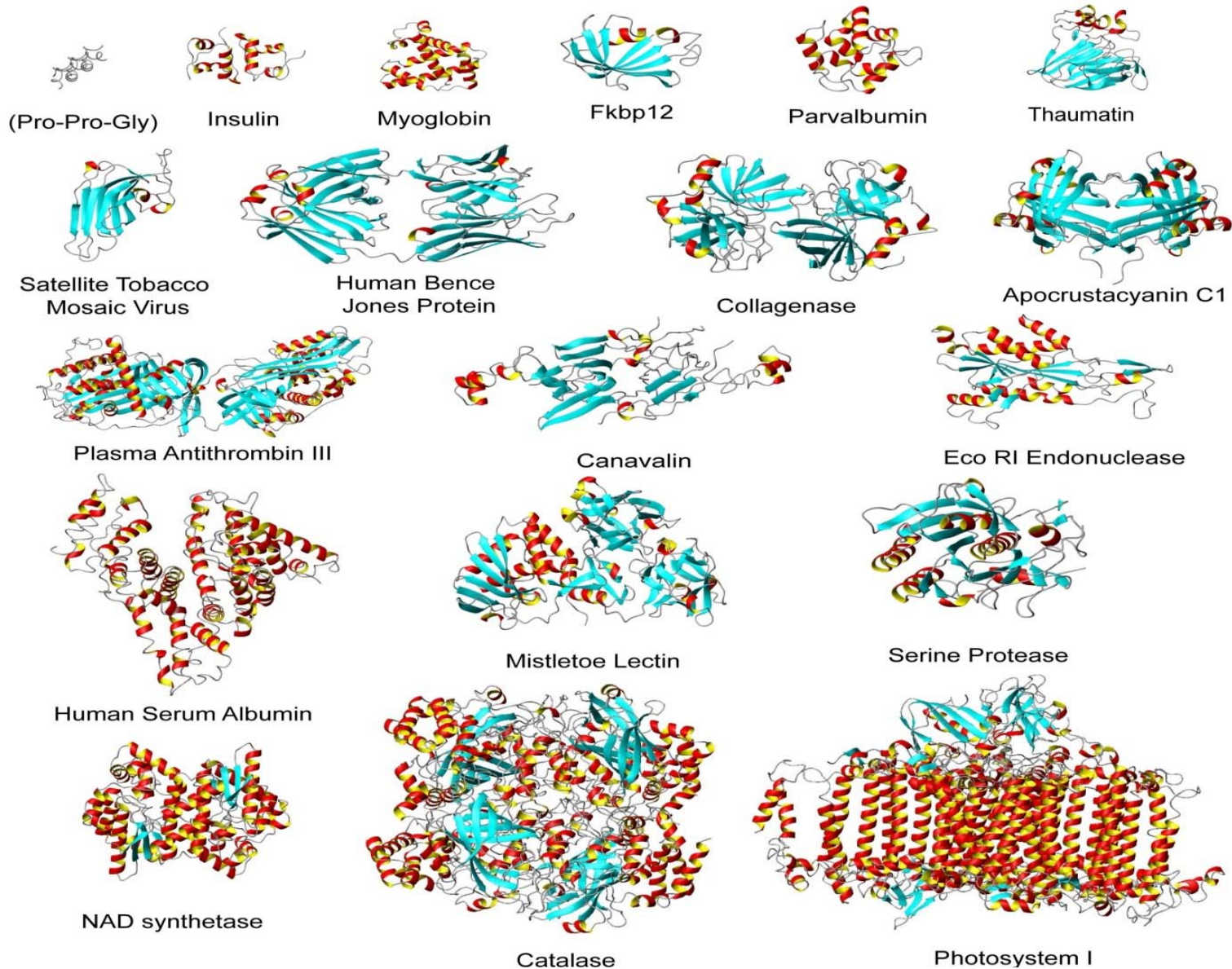

# Examples of Microgravity Grown Crystals


## *Recombinant Human Insulin*

100 X Magnification

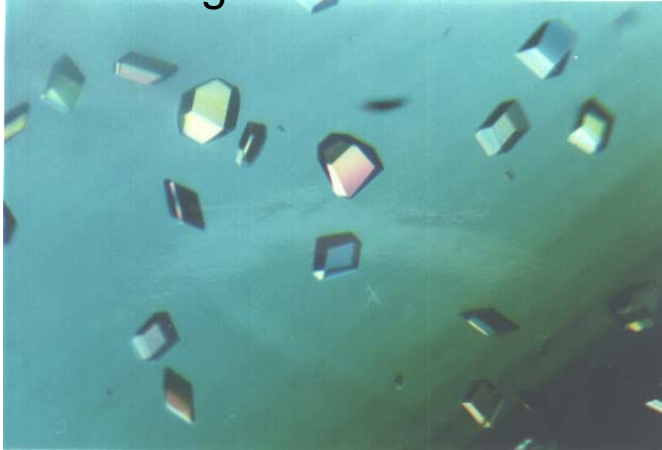


Earth-Grown Insulin


10 Million X Magnification

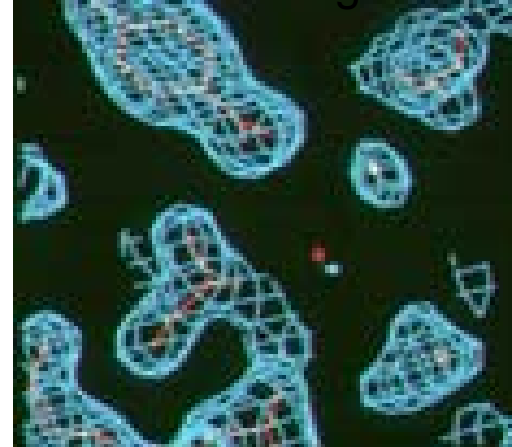

Earth-Grown Insulin

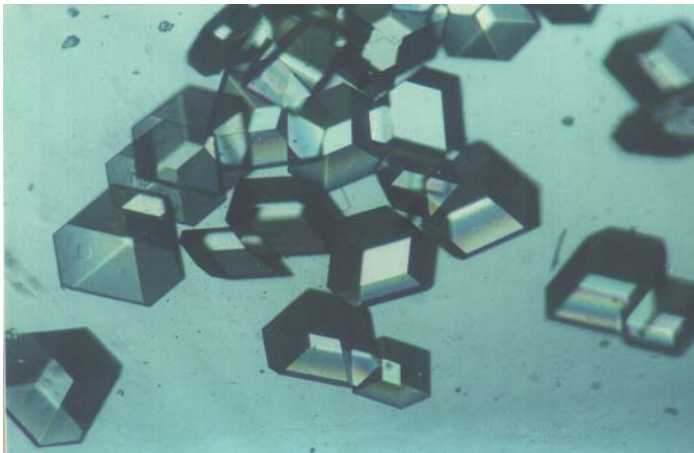

Space-Grown Insulin

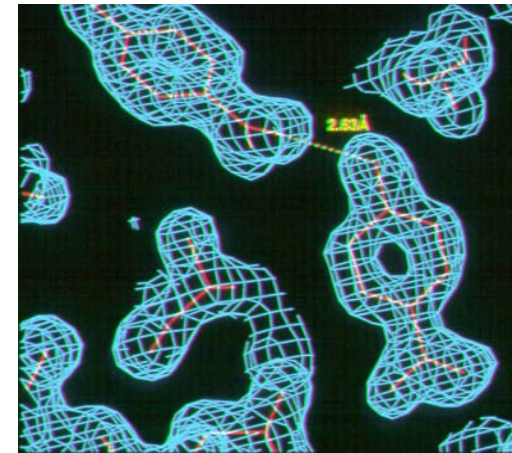

Space-Grown Insulin



# Microgravity Biomaterials Recent Selections

---



Biophysics/Biomaterials currently has 4 investigations selected from proposals received in response to a 2014 NASA Research Announcement

“Solution Convection and the Nucleation Precursors in Protein Crystallization”  
by Dr. Peter Vekilov of the University of Houston

“Growth Rate Dispersion as a Predictive Indicator for Biological Crystal Samples Where Quality Can be Improved with Microgravity Growth” by Dr. Edward Snell of the Hauptmann Woodward Medical Institute

“The Effect of Macromolecular Transport on Microgravity Protein Crystallization”  
by Dr. Lawrence Delucas of the University of Alabama at Birmingham

“Amyloid Fibril Formation in Microgravity: Distinguishing Interfacial and Flow Effects” by Dr. Amir Hirsu of Rensselaer Polytechnic Institute

Note that the Center for the Advancement of Science in Space (CASIS) is also funding ISS researchers in this area



# Microgravity Biomaterials Workshop

---



Future: To determine the level of interest in new areas of microgravity biomaterials research a Request for Information was released by NASA

- 23 responses

In addition, a workshop was held on December 3 at the recent 2014 Fall Materials Research Society Meeting and Exhibition in Boston entitled, "NASA Biological Materials, Biomaterials and Biomimetics Workshop" to further evaluate the level of interest in conducting biomaterials research on ISS.

- Approximately 50 Attendees





# Microgravity Biomaterials RFI Requested Data

---



What best describes your area of expertise?

What category is your employer?

Please describe at a high level the materials experiment or set of closely related experiments which you believe would be the highest scientific quality and the widest impact to the research community that require the use of the International Space Station.

Please describe the equipment, hardware, facilities, etc. that would be necessary to perform the experiment(s) which you described above on the ISS.

What is the application and broader impact of this research to life on earth and/or space exploration?

Which (if any) other government agencies or commercial entities would you anticipate partnering with NASA in support of performing the experiment(s) you described above?

What types of experimental data and approximately how much data, e.g. megabytes, gigabytes, etc. would be generated by the experiment(s) which you describe above?



# Microgravity Biomaterials Areas of Interest

---



Freeze casting of bioscaffolds from aqueous solutions


Electro-spinning of fibers to be used in patterned bioscaffolds

Ordered biomaterials structured from acoustic fields

Bioscaffolds constructed by sol-gel transitions



## Available Hardware: EXPRESS Rack

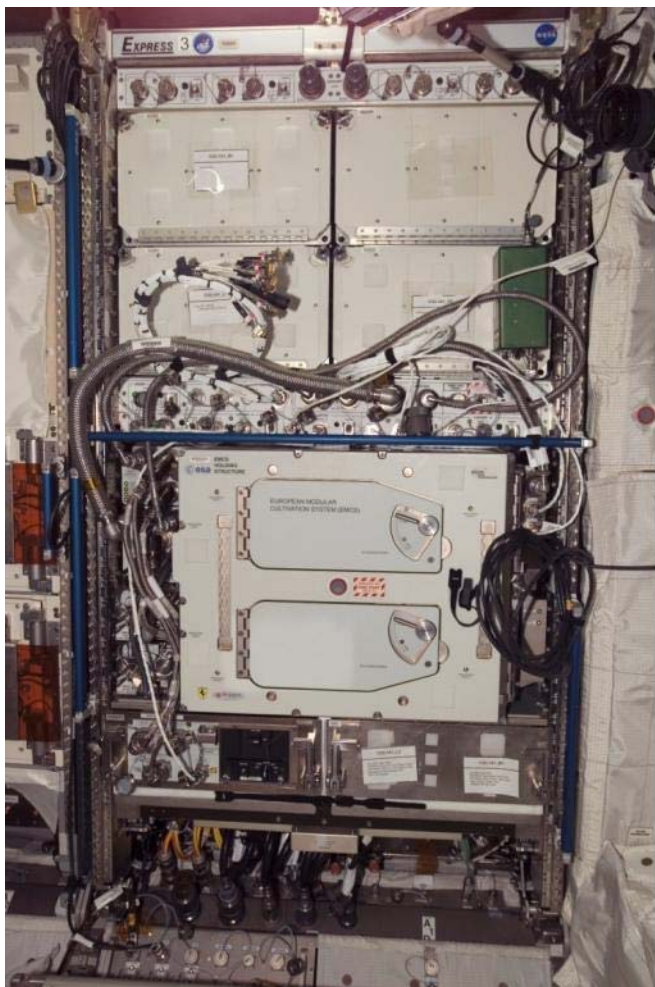

EXPRESS Rack 3 with European Modular Cultivation System operating in Columbus Module

- EXpedite the PROcessing of Experiments to Space Station (EXPRESS) rack is a multi-use facility which provides standard interfaces and resources for Middeck Locker-type and International Subrack Interface Standard (ISIS) Payloads
  - 28Vdc power
  - Ethernet, RS422 , Analog, Discrete
  - Air and Water (2 locations per rack) Cooling
  - NTSC Video
  - Vacuum Exhaust (1 location per rack)
  - Nitrogen Supply (1 location per rack)
- Eight flight racks
- Trainer Racks at JSC and MSFC to support crew and ground training
- Functional Checkout Unit (FCU) to support payload testing at MSFC





# EXPRESS Payload Resources




Resource	Amount per Payload Position	
	Locker	ISIS Drawer
<b>Structural Attachment</b>	Attachment to Rack per IDD •Mass constraint launch vehicle dependent	Attachment to Rack per ISIS Spec •64 lb within cg constraints
<b>Power</b>	5, 10, 15, or 20 Amp at 28 VDC	5, 10, 15, or 20 Amp at 28 VDC
<b>Thermal Control Air</b>	Nominal 150 W (1200 W rack maximum)	Nominal 150 W (1200 W rack maximum)
<b>Thermal Control Water</b>	500 W Heat Rejection per position (2 positions per rack)	500 W Heat Rejection per position (2 positions per rack)
<b>Data</b>	•1 - RS-422    •2 - +/- 5 Vdc Analog •1 - Ethernet    •3 - 5 Vdc Discrete (bi-dir)	•1 - RS- 422    •1 - +/- 5 Vdc Analog •1 - Ethernet    •2 - 5 Vdc Discrete (bi-dir)
<b>Video</b>	NTSC/RS 170A feed from payload source (Shared)	NTSC/RS 170A feed from payload source (Shared)
<b>Venting</b>	1 payload interface per rack (Shared)	1 payload interface per rack (Shared)
<b>Nitrogen</b>	1 payload interface per rack (Shared, 12 lbm/hr)	1 payload interface per rack (Shared, 12 lbm/hr)

Reference: EXPRESS Rack Payloads Interface Definition Document, SSP 52000-IDD-ERP



# EXPRESS Rack Drawer

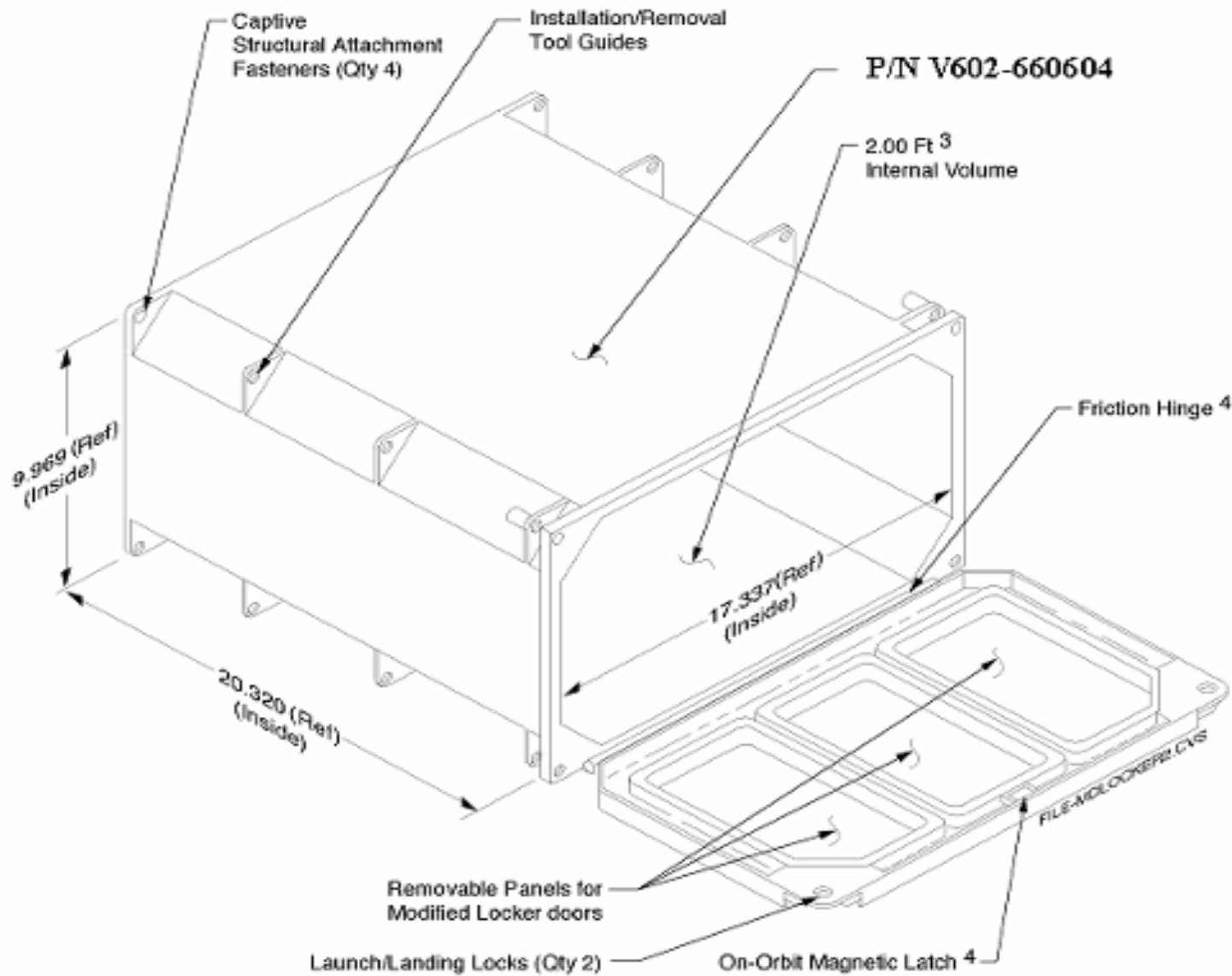



# Available Hardware: Microgravity Science Glovebox



## Removable Side Ports

16" diameter on both Left and Right sides for setting up hardware in Work Volume


## Glove Ports

Four identical glove ports are located on the left and right side loading ports and the front window

## DC Power Switching And Circuit Breakers

## Stowage Drawers

## Video System Drawer

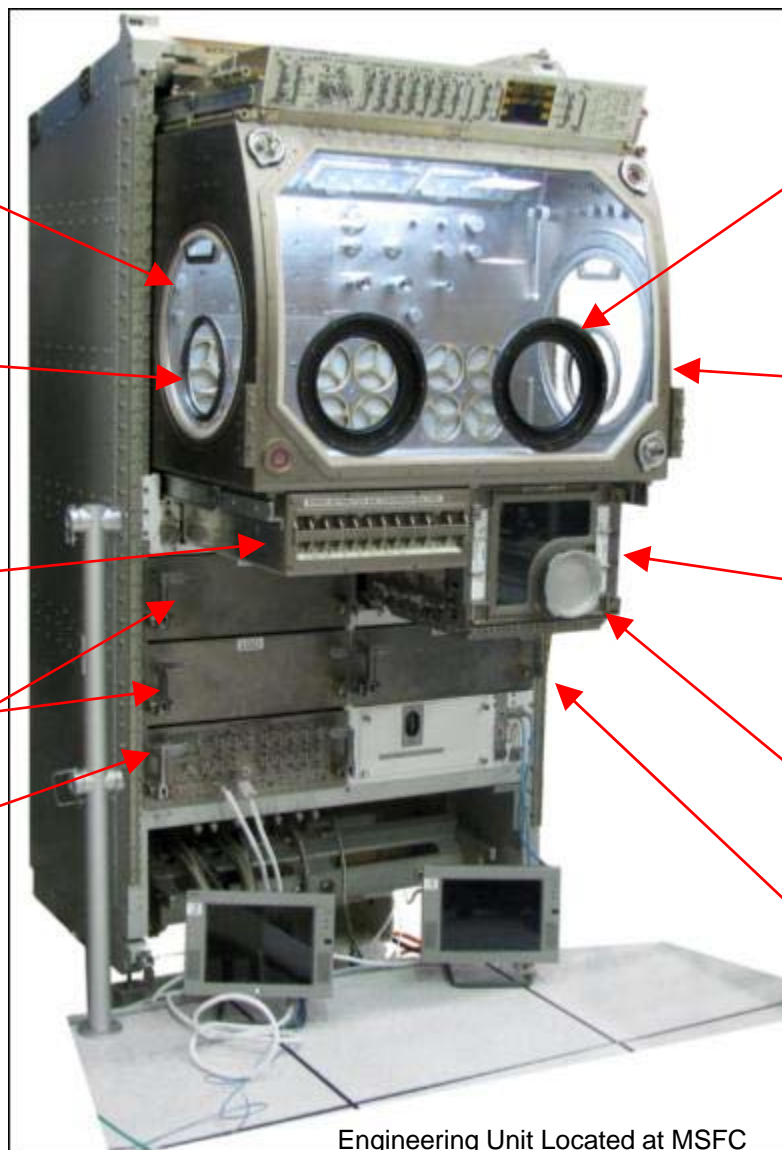

## Front Window Glove Ports

Four 6" diameter glove ports can be fitted with any of three different sized gloves or blanks

## Core Facility

Retractable Core Facility includes the Work Volume, Airlock, Power Distribution & Switching Box, and the Command and Monitoring Panel

## Airlock

Provides a "Pass Through" for hardware to enter the Work Volume without breaking Containment. The lid of the Air Lock opens up into the floor of the Work Volume

## Airlock Glove Port with Blank

A Single 4" diameter glove port can also be fitted with any of three different sized gloves or a blank

## Stowage Drawers



# Current MSG-Provided Payload Interfaces/Resources

---




- **Work Volume(WV) - Volume**
  - 0.255 m<sup>3</sup> = 255 liters
- **Work Volume - Dimensions**
  - 906mm wide x 637mm high
  - 500mm deep (at the floor)
  - 385mm deep (at the top)
- **Maximum size of single piece of equipment in WV (via side access ports)**
  - 406mm diameter
- **Payload Attachment**
  - M6 threaded fasteners in floor, ceiling, & sides
- **Power available to investigation**
  - +28V DC at useable 7 amps
  - +12V DC at useable 2 amps
  - -12V DC at useable 2 amps
  - +5V DC at useable 4 amps
  - +120V DC at useable 8.3 amps
- **Maximum heat dissipation**
  - 1000W Total
    - 800W from coldplate
    - 200W from air flow
- **General illumination**
  - 1000 lux @ 200mm above WV floor
- **Video**
  - 4 color Hitachi HV-C20 cameras
  - 2 Sony DSRV10 Digital Recorders
  - 2 Sony GV-A500 Analog 8mm Recorders
- **Data handling connections**
  - Two RS422-to-MSG for investigations
  - One MIL-BUS-1553B-to-MSG for communication via MLC
  - Ethernet LAN 1 and LAN 2 (in US LAB)
  - MSG Laptop Computer (MLC) – IBM T61P
- **Filtration**
  - 12 HEPA/charcoal/catalyst WV filters
- **1 HEPA/charcoal/catalyst Airlock filter**
- **Up to Two Levels of Containment**
  - Physical barrier of MSG structures, gloves, etc.
  - Negative pressure generated by MSG fans.
- **Other resources available**
  - Gaseous Nitrogen
  - Vacuum (VRS & VES)





# Available Hardware: Light Microscopy Module

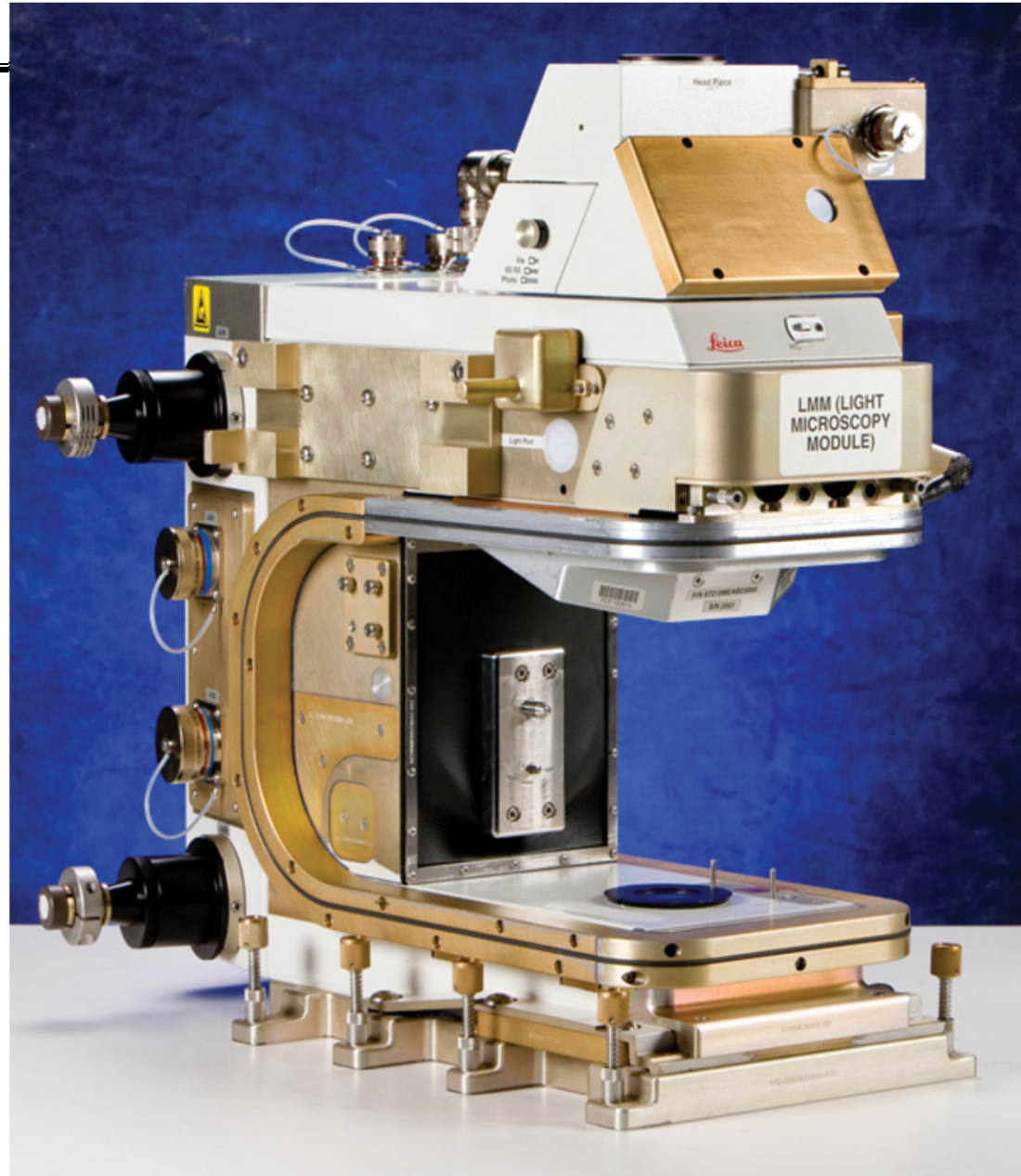



# Light Microscopy Module

---



The microscope can house many different lenses corresponding to magnifications of 2.5, 4, 10, 20, 40, 50, 63 (air), 63 and 100 oil-coupled objectives.

Present capabilities include brightfield, fluorescence, and epi-illumination microscopy.

Future planned capabilities include high-resolution color video microscopy, condenser assembly, confocal microscopy, and possibly laser tweezers.

# Available Hardware: Fluorescence Microscope


## Specifications

- Inverted microscope: Leica DMI 6000B
- Bright field, phase contrast, and fluorescence observations
- Objective lens: x5 (NA 0.12), x10 (NA 0.25), x20 (NA 0.35), x20 (NA 0.70), x40 (NA 0.50), x40 (NA 0.75)
- Monochrome CCD camera: Leica 360FX
- GFP and DsRED observations
- LED illumination (365, 470, 530, 620nm): Leica SFL7000
- Time-lapse videomicroscopy
- Leica micro-titer plate holder (No. 11531434) is set as a standard specimen holder.
- The maximum size of the specimen vessel including frame or holder is 83x127mm (micro-titer plate size)
- Objective lens, CCD camera, fluorescence filter cube, and LED illumination are exchangeable to Leica series

Launched with HTV-3 in 2012, and onboard check out has been completed.



JAXA Microscope in MSPR

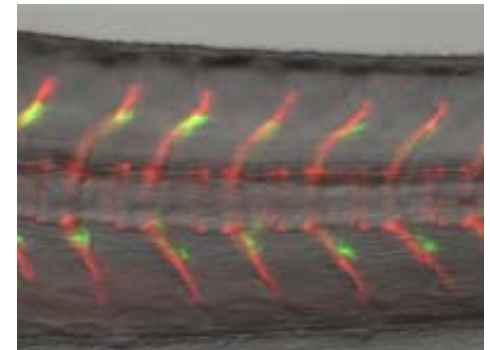

Medaka transgenic line



## Other Hardware

---



A number of refrigerators, freezers, incubators, and crystal growth hardware devices also are available.





# Microgravity Biomaterials Wrap Up

---



Next Step: Likely to release an NRA in the next year.

Finally –

Comments, questions, or inputs on the program are welcome

Point of Contact: James Patton Downey

Organization: EM31/MSFC/NASA

Phone 256-544-6432

Email: [patton.downey@nasa.gov](mailto:patton.downey@nasa.gov)